



rubica



CATALOG OF DEVICES

Rubica devices:

Smart PM	3
<i>Modern device for micropigmentation</i>	
CapilMed Electrocoagulator	4
<i>Device for sealing broken vessels and removing unwanted hair</i>	
Sonic Peel	5
<i>Cavitation spatula</i>	
Mezolift Massager	6
<i>Massager 5 in 1</i>	
Photon LED Skin Lamp	7
<i>Face and body photon therapy lamp</i>	
Multi Massager 4D	8
<i>Multifunctional treatment massager</i>	
IQ Skin Massager	9
<i>Smart face massager</i>	
Compact H3 Pen	10
<i>Device for microneedle mesotherapy with simultaneous cocktail administration</i>	
Electroporation Pen	11
<i>Device combining microneedling with microcurrents and LED light therapy</i>	
Darsonval Rubica	12
<i>Darsonval device with 4 electrodes</i>	
Vacufrax RF	13
<i>Vacuum microneedle radiofrequency device</i>	
Mesopen	14
<i>Device for microneedle mesotherapy treatments</i>	
Eye & Lip Beauty Massager	15
<i>Eye and lip area massage device</i>	

CODE: 7046
WARRANTY: 3 YEARS

Smart PM

Next-generation device for permanent make-up, micro-needle mesotherapy and micropuncture

Advantages of the device:

- Polish production
- modern design
- ergonomic aluminum handle
- intuitive control panel made of tempered glass
- smooth speed regulation
- precise control of needle ejection
- quiet, precise operation
- software in 5 languages: Polish, English, German, Spanish, Russian

Micropigmentation modules:

- traditional needlesh
- acupunture needles
- powder method needles
- powder method modules



Treatment options:

- permanent makeup
- medical micropigmentation
- microneedling mesotherapy
- micropuncture



TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Size	129 x 73 x 95 mm
Station weight	345 g
Handle and cable weight	101 g
Handle stand weight	21 g
Power	24 W
Power supply	100 - 240 AC / 12DC
Needle ejection speed	50 - 150 Hz adjustable scale

The set includes:

- main unit
- handle
- protective case with handle and strap
- power supply
- foot switch
- handle stand
- manual and warranty card
- extension cable for the handle



CODE: 7048
WARRANTY: 2 YEARS

CapilMed Electrocoagulator

A modern, 3-function device with 100% effectiveness in closing broken vessels and removing unwanted hair

Advantages of the device:

- Polish manufacturing
- modern design
- intuitive touch panel
- precise selection of parameters
- stainless steel needles
- Software in 5 languages: Polish, English, German, Spanish, Russian

Methods of work using the device:

- thermolysis
- electrolysis
- blend method

Treatment options:

- hair removal
- fibroma removal
- closing of dilated blood vessels




The set includes:

- active electrode handle with needle holder
- passive electrode
- foot switch
- 3-needle set
- user's manual with the warranty card

AVAILABLE CARTRIDGES

Ballet K2 electrocoagulation needle 1 pc.	70481
Ballet K3 electrocoagulation needle 1 pc.	70482
Ballet K4 electrocoagulation needle 1 pc.	70483

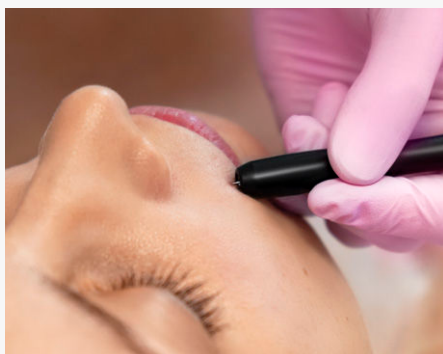

Suggested treatment:

SEALING OF VESSELS AND REDUCTION OF ERYTHEMA

treatment that eliminates the causes and symptoms of rosacea

- Step 1.** Remove make-up with **Micellar Vit. U Lotion**.
- Step 2.** Perform enzymatic peeling with **Sensi Enzymatic Peel**.
- Step 3.** Tone the skin with **Vit. U Tonic**.
- Step 4.** Massage in **Vit. U Stop Couperose Serum**.
- Step 5.** Apply **Hydroalgae Vit. U Mask**.
- Step 6.** Apply **Vit. U Anti-Couperose Mask** to the neck and cleavage.
- Step 7.** Apply **Vit. U Cream**.

CODE: 7049
WARRANTY: 12 MONTHS

Sonic Peel

Ultrasound spatula

Methods of work using the device:

- **cavitation peeling** - gently removes dead skin cells, leaving the skin smooth and radiant
- **cavitation sonophoresis** - helps improve the absorption of active cosmetic ingredients
- **iontophoresis (+) and (-)** - regulates electrical charges to improve skin condition



Suggested treatment:

BAKUCHIOL CLEANSING

deep cleansing treatment and normalization of sebaceous glands

Step 1. Remove make up with **Max Matt Gel**.

Step 2. Perform enzymatic peeling with **AHA Max Forte Peeling**.

Step 3. Tone the skin using **Max Matt Tonic**.

Step 4. Soften the skin with **Max Derma Soft Mask**.

Step 5. Perform cavitation peeling using **Sonic Peel**.

Step 6. Perform darsonvalization with **IQ Digital Darsonval**.

Step 7. Apply **Max Sebum Control Serum** and perform cavitation sonophoresis.

Step 8. Apply **Max Balance Mask**.

Step 9. Apply **Max Sebum Control Cream**.



TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Package dimensions	205 x 75 x 40 mm
Product weight	78g
Battery voltage	3,7 V
Rated power	1,5 W
Battery capacity	500 mAh
Blade oscillation frequency	24 kHz $\pm$ 1 kHz
Run time	4 minutes
No cable operation capability	yes
Maximum operating time on a single charge	2 hours

The set includes:

- device,
- charging cable,
- instruction manuals



CODE: 7086
WARRANTY: 12 MONTHS

Mezolift Massager

Needle-free mesotherapy,
5-in-1 massager, intended for home use

The device is used to perform cleansing, moisturizing, regenerating and anti-aging treatments. The multifunctional device provides non-invasive skin care for the face, neck and décolleté and better absorption of cosmetics.

Methods of work using the device:

- needle-free mesotherapy
- EMS microcurrents
- phototherapy
- HOT-COLD treatment
- vibration micromassage

Treatment options:

- cleansing treatments
- moisturizing treatments
- regenerating treatments
- anti-aging treatments



TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Size	175 x 55 x 67 mm
Working programs	5
Power	8,5 W
Battery capacity	1000 mAh
Charging voltage	500 mA
Charging time	± 2,5 hours

The set includes:

- main unit
- USB cable
- stand
- manual
- box

Suggested treatment:

EXOSS PRO-AGE LINE WITH MEZOLIFT MASSAGER

exosomal-senolytic therapy in combination with
Mezolift Massager

- Step 1.** Remove make-up with **Caviar Micellar Gel**.
- Step 2.** Massage the neck with the **Mezolift Massager**.
- Step 3.** Perform peeling using **Caviar Face Peeling**.
- Step 4.** Tone the skin using **Caviar Micellar Tonic**.
- Step 5.** Apply **ExoLytic Fusion Serum** and massage it in using the device.
- Step 6.** Apply **Mousse-Oil Mask**.
- Step 7.** Apply **ExoLytic Fusion Cream**.

CODE: 7087
WARRANTY: 12 MONTHS

Photon LED Skin Lamp

Photon therapy lamp for face and body

Advantages of the device:

- 6 colors and a variable cycle
- 3 intensity levels of each color
- timer 0-60 minutes, with automatic switch off (suggested irradiation time is 15 min)
- temperature adjustment 30-70°C
- protective goggles included
- 100 LED points guarantee even illumination

Treatment options depending on the light colour:

- **violet** - supports the process of wound healing
- **blue** - soothes irritation and inflammation, indication for acne-prone skin
- **turquoise** - reduces inflammation of capillaries, relaxes the skin
- **green** - reduces seborrhea, improves color and reduces the tendency to discoloration
- **yellow** - soothes redness, smooths and lifts the skin
- **red** - stimulates microcirculation, regeneration and collagen stimulation

Suggested treatment:

ULTRASOUND PHOTON THERAPY

exosomal-senolytic needle-free mesotherapy and light therapy for mature skin

Step 1. Remove make-up with **Caviar Micellar Gel**.

Step 2. Perform peeling using **Caviar Face Peeling**.

Step 3. Tone the skin with **Caviar Micellar Tonic**.

Step 4. Apply **ExoLytic Fusion Serum**.

Step 5. Apply **Sensi Peptide Sono Gel**.

Step 6. Perform needle-free mesotherapy.

Step 7. Perform photon therapy with **Photon LED Skin Lamp**.

Step 8. Perform nourishing and stimulating massage.

Step 9. Apply **ExoLytic Fusion Cream**.



TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Model	BM-NFY7C
Dimensions	435 x 265 x 329 mm
Weight	4.5 kg ± 10%
Power supply	DC24V4A50HZ/60HZ
Power	96 W
Lifespan of the photon lamp	10,000 hours
LED Color Light	6 colors

The set includes:

- main unit
- power cable
- instruction manual
- protective goggles
- box



CODE: 7088
WARRANTY: 12 MONTHS

Multi Massager 4D

Multifunctional treatment massager

Advantages of the device:

- 4 vibrating-massaging rollers, rotating 360 degrees
- perfectly adapts to the massaged area
- stimulates skin, subcutaneous tissue and muscles
- cordless
- with a convenient, aesthetic charging base
- combines massage with the introduction of active ingredients

Methods of work using the device:

- 4D massage
- EMS
- RF
- photon therapy - 4 colors
- heat/cold



Suggested treatment:

MULTI-LEVEL MOISTURIZATION WITH 4D MASSAGE

hyaluronic treatment combined with Multi Massager 4D

- Step 1.** Remove make-up with **Hydro Intense Emulsion**.
- Step 2.** Perform neck massage with **Multi Massager 4D**.
- Step 3.** Perform peeling with **Hydro Intense Peel**.
- Step 4.** Tone the skin using **Hydro Intense Tonic**.
- Step 5.** Apply **Hydro Deep Intense Serum** and massage it in using the massager.
- Step 6.** Apply **Hydro Intense Gel Mask** to the face, neck and cleavage.
- Step 7.** Apply **Hydro Life Intense Cream**.

CODE: 7089

WARRANTY: 12 MONTHS

IQ Skin Massager

Smart face massager

Advantages of the device:

- vibration massage, EMS and LED light therapy
- improves muscle tone and reduces wrinkles
- stimulates collagen production and blood circulation
- enhances skin nourishment and firmness
- designed for safe and easy home use
- suitable for all skin types

Methods of work using the device:

- **clean mode** - deeply cleanses and narrows pores, improves skin tone. This mode uses **green** light.
- **EMS mode (microcurrent)** - introduces the active ingredients of products deep into the skin. It is recommended to use anti-aging, moisturizing serums. This mode uses **red** light.
- **hot mode (heat)** - massage using pleasant heat improves circulation, removes tension, firms the skin. This mode uses **blue** light
- **micro vibrations (vibration massage)** used in each of the three modes open the skin channels, which improves the absorption of active ingredients.

Treatment options depending on the light colour:

- LEDs emit light of different wavelengths with no thermal effect:
- **red light** - there is a dilation of blood vessels, improves microcirculation, stimulates collagen synthesis, facilitates the exchange of nutrients
- **green light** - improves blood circulation, reduces sebum secretion, accelerates lymph flow
- **blue light** - has antibacterial, bactericidal effects, soothes inflammation, works superficially at the epidermal level



Suggested treatment:

MENOSHE TREATMENT

reduces the symptoms of skin aging during menopause

- Step 1.** Remove make-up with **Hydro Intense Emulsion**.
- Step 2.** Perform peeling with **Hydro Intense Peel**.
- Step 3.** Tone the skin with **Hydro Intense Tonic**.
- Step 4.** Apply **MenoShe Serum**.
- Step 5.** Perform massage with **IQ Skin Massager**.
- Step 6.** Apply **MenoShe Mask** to the neck and cleavage.
- Step 7.** Apply **Algaplast MenoShe Mask**.
- Step 8.** Apply **MenoShe Cream**.

CODE: 7093
WARRANTY: 12 MONTHS

Compact H3 Pen


Device for microneedle mesotherapy with simultaneous cocktail administration

Advantages of the device:

- cartridges integrated with a cocktail dispenser
- LED display with 6-level adjustment
- needle extension adjustment from 0.25 -1.0 mm
- cordless
- light and mobile
- treatments for the face, neck, cleavage, hands and body

Treatment options

- wrinkle reduction
- skin thickening and firming
- color uniforming
- reducing the visibility of scars and discolorations
- improving the contour of the face

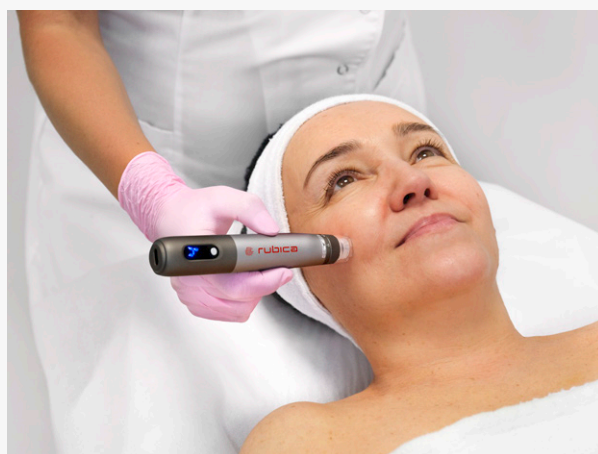

TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Dimensions	169 x 123 x 56 mm
Weight	390 g
Input power	100 - 240 V
Output power	5 V

The set includes:

- device
- power cable
- user's manual
- adapter
- 3 cartridges: H12, H24, nanoHR

AVAILABLE CARTRIDGES

Cartridge 12 needle Compact H3 Pen 1 pc.	70931
Cartridge 24 needle Compact H3 Pen 1 pc.	70932
Cartridge nano needle Compact H3 Pen 1 pc.	70933
Adjusting ring for Compact H3 Pen 1 pc.	70934

Suggested treatment:

COMPACT NANOMESOTHERAPY H3 PEN

multi-level moisturization combined with peptide nano-peeling

- Step 1.** Remove make-up with **Hydro Intense Emulsion**.
- Step 2.** Cleanse the skin with **Pre Peel Cleanser**.
- Step 3.** Fill in the cartridge
- Step 4.** Perform peptide nanomesotherapy.
- Step 5.** Neutralize the skin with **Neutralizer Lotion**.
- Step 6.** Tone the skin using **Hydro Intense Tonic**.
- Step 7.** Apply **Hydro Deep Intense Serum**.
- Step 8.** Apply **Hyaluron Crystal Collagen Mask**.
- Step 9.** Apply **Hydro Intense Gel Mask** to the neck and cleavage.
- Step 10.** Apply **Hydro Ever Cream**.

CODE: 7094
WARRANTY: 12 MONTHS

Electroporation Pen


The device combines microneedling with microcurrents (EMS) and LED light therapy

Advantages of the device:

- wireless
- mobile
- digital display
- 5-step speed control
- 0.1-1.5 mm needle extension settings
- the induction base stops the working device and resumes operation after it is lifted

Methods of work using the device:

- microneedle mesoporation
- microneedle mesotherapy
- LED: red and blue

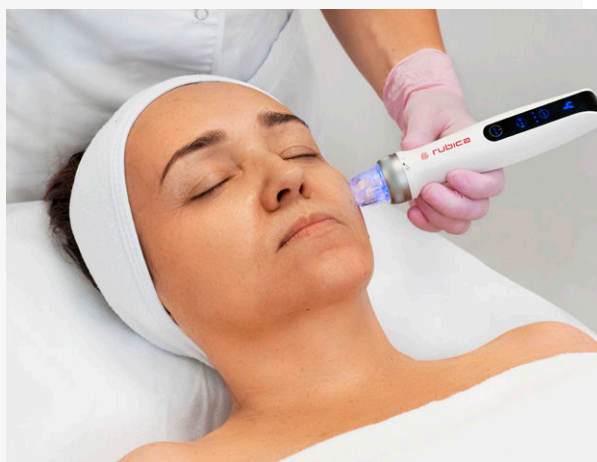

TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Dimensions	193 x 146 x 70 mm
Weight	670 g
Input power	100 - 240 V
Output power	5 V

The set includes:

- device
- induction base
- adapter
- power cable
- 3 cartridges: nano, 14-needles and 25-needles
- user's manual

AVAILABLE CARTRIDGES

Electroporation Pen 14 needle module 1 pc.	70941
Electroporation Pen 25 needle module 1 pc.	70942
Electroporation Pen nano needle module 1 pc.	70943

Suggested treatment:

EXOSS PRO-AGE TREATMENT

exosomal-senolytic microneedle mesoporation of mature skin

- Step 1.** Remove make-up with **Sensi Peptide Foam**.
- Step 2.** Perform peeling with **Sensi Enzymatic Peel**.
- Step 3.** Tone the skin with **Sensi Peptide Tonic**.
- Step 4.** Apply **ExoLytic Fusion Serum**.
- Step 5.** Perform microneedle mesoporation.
- Step 6.** Apply **Konjac Fiber Mask**.
- Step 7.** Apply **ExoLytic Fusion Cream**.

CODE: 72010
WARRANTY: 12 MONTHS

Darsonval Rubica


Set with Darsonval device
and 4 electrodes

The kit comprises a device (main unit) and 4 glass wands with the following attachments:

- tongue electrode
- curved tip electrode
- mushroom electrode
- comb electrode with Tungsten

Treatment options:

- inhibits sebum secretion
- stimulates exfoliation of dead skin
- firms and smooths
- stimulates skin regeneration and renewal processes
- reduces wrinkles
- improves blood supply to the skin
- reduces dilated pores
- has disinfectant properties
- removes acne lesions
- inhibits the progression of baldness
- facilitates herpes healing

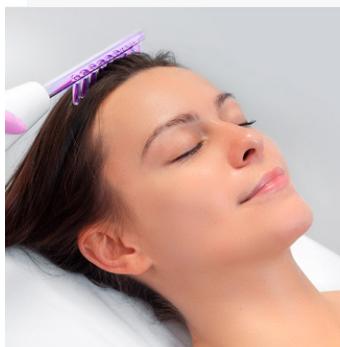

The set includes:

- main unit
- power cable
- 4 glass wands (tongue electrode, curved tip electrode, mushroom electrode, comb electrode with Tungsten)
- user's manual

Suggested treatment:

BAKUCHIOL CLEANSING

treatment for deep cleansing and normalization of the skin

- Step 1.** Remove make-up with **Max Matt Gel**.
- Step 2.** Perform peeling with **AHA Max Forte Peeling**.
- Step 3.** Tone the skin using **Max Matt Tonic**.
- Step 4.** Soften the skin with **Max Derma Soft Mask**.
- Step 5.** Perform cavitation peeling and manual cleansing.
- Step 6.** Perform darsonvalization with **Darsonval Rubica**.
- Step 7.** Apply **Max Sebum Control Serum** and perform cavitation sonophoresis.
- Step 8.** Apply **Max Balance Mask** for 20 minutes.
- Step 9.** Apply **Max Sebum Control Cream**.

CODE: 7095
WARRANTY: 12 MONTHS

Vacufrax RF


Micro-needle radiofrequency device that combines, micro-puncture with vacuum and radio waves

Advantages of the device:

- menu in 8 languages: Polish, English, German, Spanish, French, Chinese, Italian, Russian
- intuitive control panel
- defined parameters
- treatment comfort
- access to hard-to-reach areas
- face and body treatments
- spectacular treatment results

Unique features:

- option to choose modules for the treatment area: eye area, wrinkles, furrows, face, and body (3 modules with different numbers of needles and treatment areas)
- additional nano module (needle-free) – allows working on particularly sensitive and vascular skin
- vacuum function – enhances client's comfort and treatment efficacy
- additional cooling head – permits faster restoration of the skin's thermal balance
- modern and legible interface with division into treatment areas

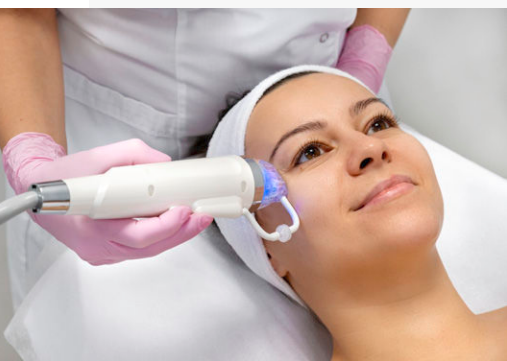

TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Maximum depth	38 cm
Approximate height	33 cm
Device width	17 cm
Width with handles	38 cm

The set includes:

- main device
- treatment head
- cooling head
- 4 modules in the device package
- user manual with warranty card

AVAILABLE CARTRIDGES

Vacufrax RF 64-needle module 1 pc.	70951
Vacufrax RF 25-needle module 1 pc.	70952
Vacufrax RF 10-needle module 1 pc.	70953
Vacufrax RF nano module 1 pc.	70954
Vacufrax RF vacuum hose 1 pc.	70955

Suggested treatment:

VACUUM MICRONEEDLE RADIOFREQUENCY

rejuvenation treatment with Vacufrax RF

- Step 1.** Remove make-up with **Micellar Vit. U Lotion**.
- Step 2.** Tone the skin with **Vit. U Tonic**.
- Step 3.** Disinfect the skin with an appropriate product.
- Step 4.** Perform microneedle radiofrequency with **Vacufrax RF**.
- Step 5.** Use the cooling head.
- Step 6.** Apply **Snail Mucin Cream**.


CODE: 7581
WARRANTY: 12 MONTHS

Mesopen

Device for microneedle mesotherapy treatments

Advantages of the device:

- stage regulation of speed levels
- precise adjustment of needle length (0.25-2 mm)
- cordless, portable
- aluminum housing
- lightweight (weight: 56 g)

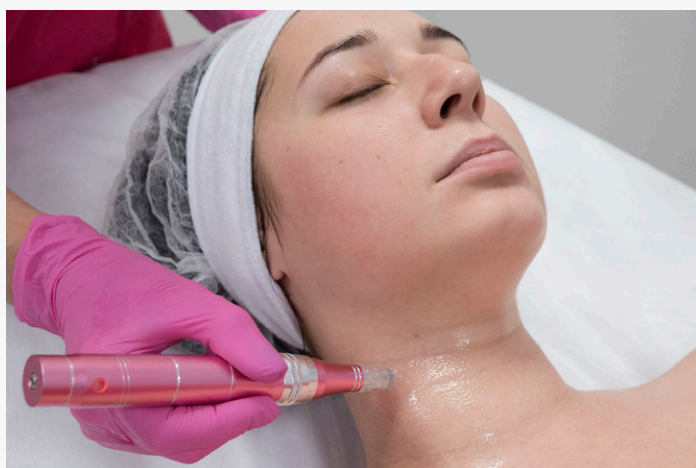

TECHNICAL PARAMETERS AND FEATURES OF THE DEVICE

Weight	56 g
Needle length adjustment	0,25 - 2 mm
Housing	aluminium

AVAILABLE CARTRIDGES

Cartridge nano for Mesopen 1 pc.	75811
Cartridge 24-needles for Mesopen 1 pc.	75812
Cartridge 12-needles for Mesopen 1 pc.	75813
Cartridge 3D for Mesopen 1 pc.	75818
Cartridge 36-needles for Mesopen 1 pc.	75819

Suggested treatment:

WHITENING NANO EXFOLIATION

treatment combining nanoneedle mesotherapy with 40% mandelic acid

- Step 1.** Remove make-up using **Sensi Peptide Foam**.
- Step 2.** Cleanse the skin with **Pre Peel Cleanser**.
- Step 3.** Perform nanoneedle mesotherapy with **Mesopen**.
- Step 4.** Perform exfoliation **MandeliC Hybrid Peel**.
- Step 5.** Neutralize the skin with **Neutralizer Lotion**.
- Step 6.** Apply **Sensi Peptide Mask**.
- Step 7.** Apply **Sensi Peptide Cream**.

CODE: 7096
WARRANTY: 12 MONTHS

Eye & Lip Beauty Massager

Eye and lip massage device combining high-frequency micro-vibrations and phototherapy

Advantages of the device:

- wireless, lightweight, portable
- LED: red and blue
- stimulates collagen and elastin fibres production
- enhances the absorption of cosmetics



The set includes:

- device
- power cable
- user manual

Suggested treatment:

NEUROPEPTIDE REJUVENATION OF THE EYE AREA

a treatment that reduces facial wrinkles, dark circles and swelling

Step 1. Remove make-up using **2-Phase Eye Cleanser**.

Step 2. Apply **Sensi Enzymatic Peel** to the skin around the eye, leave for 10 minutes, then massage with wet hands and rinse off.

Step 3. Tone the skin with **Vit. U Tonic**.

Step 4. Apply **Eye Microneedle Argilift Cocktail**.

Step 5. Introduce the serum using **Eye & Lip Beauty Massager**.

Step 6. Apply **Tranexamic EGF Eye Mask**.

Step 7. Apply **Argilift Eye Cream**.



Rubica devices represent the perfect blend of advanced technology and user-focused design, created to meet the high standards of today's beauty professionals. Each device is engineered to deliver precise, effective, and consistent results, helping you enhance your treatments and exceed your clients' expectations. With Rubica, you invest not only in innovative solutions but also in reliability and ease of use, ensuring smooth operation and long-term satisfaction. Elevate your professional practice and stay at the forefront of the industry by choosing Rubica – where quality meets performance.

CONTACT DETAILS



+48 512 484 494
+48 782 720 105
+48 512 044 473



export@clarena.pl



www.clarena.pro/b2b_en



Clarena
ul. Clareny 2
51-361 Wrocław-Wilczyce
POLAND

Popławska Group Spółka Jawna

51-361 Wrocław, Wilczyce; ul. Clareny 2

NIP 896-139-08-45, REGON 020104340

SR dla Wrocławia Fabrycznej, IX Wydział Gospodarczy,

Nr KRS 0000891508

Offer only for customers running their own business.

Product photos presented in the offer are illustrative. Products are intended for professional use, excluding products marked for resale. The presented price offer is for information purposes only and does not constitute a commercial offer within the meaning of Art. 66 of the Civil Code. The typesetting was closed on May 29, 2023, so data created after that date could not be included in the above content. Despite our best efforts in proofreading, we also reserve the right to change the data contained in this publication. All prices presented in this catalogue include VAT.